


NOVATRAK

- Innovative filer media completely made of polypropylene nanofibers
- High dirt holding capacity
- Sanitizable and sterilizable
- Thermowelded construction
- EC listed materials for food contact
- FDA listed material per CFR21

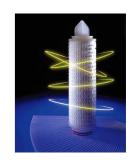

NOVATRAK cartridge represents a further milestone for Bea Technologies research and development of nanofibers for critical applications that require absolute filtration grade and reduction of the bacterial load.

NOVATRAK filter element is manufactured with an innovative filter media made of a matrix of very fine polypropylene fibers (nanofibers); they allow to achieve constant and reproducible performance, that is not influenced by the quantity of the contaminant in the liquid to be filtered, and to retain various microorganisms and bacteria as per the filtration rating.

The chemical compatibility of the materials and the technology adopted in the manufacturing process allow an effective regeneration process both by backwashing and by chemical agents.

Typical applications are in food & beverage and treatment of liquids that require particulate filtration, such as in the pharmaceutical industry.

MATERIALS OF CONSTRUCTION

Filter media	Matrix of nanofibers of polypropylene	
Upstream supports	polypropylene	
Downstream supports	polypropylene	
Internal core	polypropylene	
External cage	polypropylene	
End caps / Adapters	polypropylene	

EC DIRECTIVE FOOD SAFETY


Novatrak filter elements meet (EU) regulation 10/2011 and its subsequent amendments and regulations (EC) 1935/2004 and 1895/2005.

RECOMMENDED OPERATING CONDITIONS


- max. continuous temperature	85 °C	
- max. cumulative time of steam sterilization	40 hours at 125 °C (with cycles of 30 minutes)	
- sanitization with hot water	90 °C max	
- sanitization with chemicals	can be sanitized by standard chemical agents	
- regenerability 3% NaOH solution up to 85 °C		
- max. differential pressure	5,0 bar at 25 °C	
- recommended change out differential pressure	2,0 bar at 25 °C	

CODE	FILTRATION RATING IN LIQUIDS	BACTERIAL RETENTION FOR 10" CARTRIDGE	MAX FLOW RATE FOR 10" CARTRIDGE I/h
NVB	0,2 μm	$ \begin{array}{lll} \hbox{Oenococcus oeni} & \geq 10^6 \\ \hbox{Enterobacteriacee} & \geq 10^6 \\ \hbox{Saccharomyces} & \geq 10^{12} \\ \hbox{Lactobacillus brevi} & \geq 10^6 \\ \end{array} $	400
NVC	0,45 μm	Oenococcus oeni $\geq 10^5$ Enterobacteriacee $\geq 10^5$ Saccharomyces $\geq 10^{11}$	600
NVD	1,0 µm	Oenococcus oeni $\geq 10^3$ Enterobacteriacee $\geq 10^3$ Saccharomyces $\geq 10^{10}$	800

WATER FLOW RATE FOR 10" CARTRIDGE

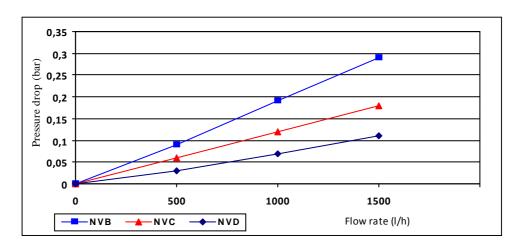

NOVATRAK ORDERING INFORMATION

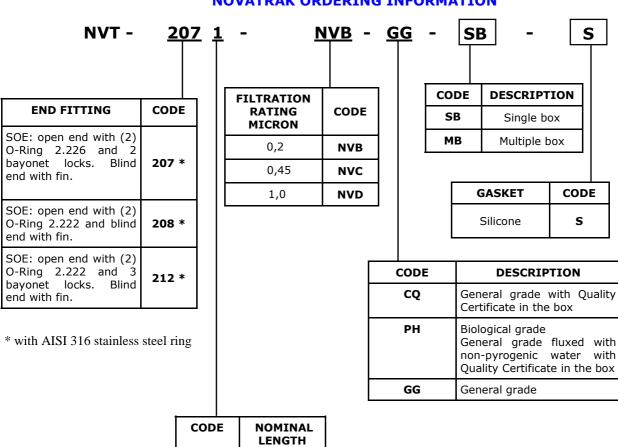

Data contained in this bulletin are informative and subject to change without notice. User is responsible for determining whether the product is fit for particular purpose and suitable for User's method of application.

10"

20"

30"

40"

1 2

3

4

