

# QUALIKAP QKP-P

## PP filter media capsules



- High filtering surface
- High mechanical resistance and dirt holding capacity
- Sanitizable and autoclave sterilizable
- EC-listed materials as per food contact
- FDA-listed materials as per CFR21
- Extractables as per USP for plastic materials
- Validation Guide available
- Low extractables



QUALIKAP QKP-P are capsules entirely manufactured using polypropylene to generate a high compatibility for fluid and aqueous solution filtration.

The pleated different filtration layers provide enhanced dirt holding capacity, together with high flow rate and a wide compatibility with the different acid/basic solutions and sanitizing agents.

The QUALIKAP QKP-P are available with filtration rating from 0,6 to 20 micron.

Manufacturing is completed in a controlled environment; each filter is subject to specific tests before delivery.

### MATERIALS OF CONSTRUCTION

<b>Filter media</b>	polypropylene
<b>Upstream drainage layers</b>	polypropylene
<b>Downstream drainage layers</b>	polypropylene
<b>Core and Cage</b>	polypropylene
<b>Terminals</b>	polypropylene
<b>Shell</b>	polypropylene

### APPLICATIONS

Ultra-pure water, alcoholic solutions, acid-base and buffer solutions, vaccines, physiological solutions, biotechnological products, ophthalmic liquids, laboratory batch purification.

### FOOD-SAFETY

QUALIKAP-QKP-P capsules meet (EU) regulation 10/2011 and its subsequent amendments and regulations (EC) 1935/2004 and 1895/2005.

### BIO-SAFETY AND EXTRACTABLES

- Materials are compliant to USP-VI CLASS toxicological requirements and USP-Plastic Materials chemical and physical requirements.
- Capsule filters meet USP "Water for injection" requirements for endotoxin particle release; the bacterial endotoxin are determined using LAL Test.
- Extractable NVR (gravimetric) after autoclave  $\leq 2$  mg.
- TOC and conductivity according USP "Purified water" and "Water for Injection" requirements.





### RECOMMENDED OPERATING CONDITIONS

- max. continuous temperature	40 °C
- max. cumulative time of steam sterilization (no in line sterilization allowed)	40 hours at 125 °C or 20 hours at 135 °C (30' cycle)
- max. pressure	5 bar at 40°C (liquids) - 3 bar at 40°C (gas)
- chemical sanitization	compatible with a wide range of sanitizers
- max. differential pressure @ forward flow	4,0 bar at 40 °C - 2,0 bar at 40°C (gas)
- max. differential pressure @ reverse flow	1,0 bar at 40 °C
- changeover recommended pressure drop	2,0 bar at 40 °C


CODE	FILTRATION RATING *	MAX SUGGESTED WATER FLOW RATE(I/h) SIZE 2	MAX SUGGESTED WATER FLOW RATE(I/h) SIZE 4
<b>KD</b>	0,6 µm	200	300
<b>KF</b>	1,2 µm	400	600
<b>KG</b>	2,5 µm	1000	1500
<b>KL</b>	4,5 µm	1500	2000
<b>KR</b>	6,5 µm	1500	2000
<b>KS</b>	10,0 µm	1500	2000
<b>KT</b>	20,0 µm	1500	2000
<b>KV</b>	40,0 µm	1500	2000
<b>KZ</b>	60,0 µm	1500	2000

\* Liquids and wet gases

## WATER FLOW RATE CHARTS

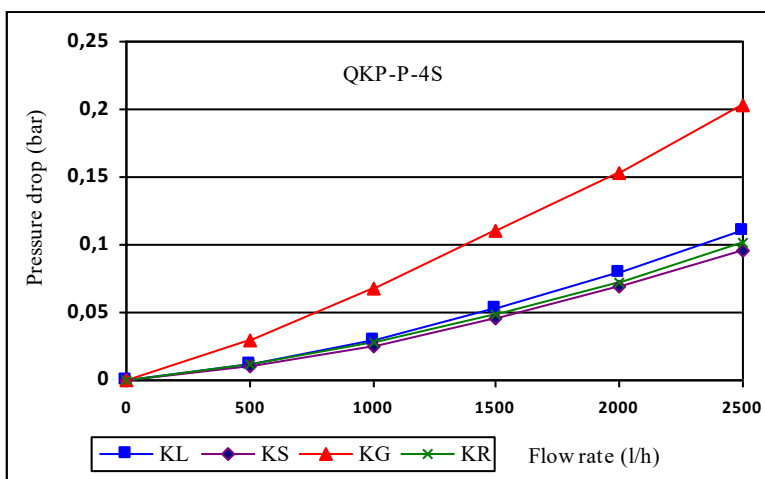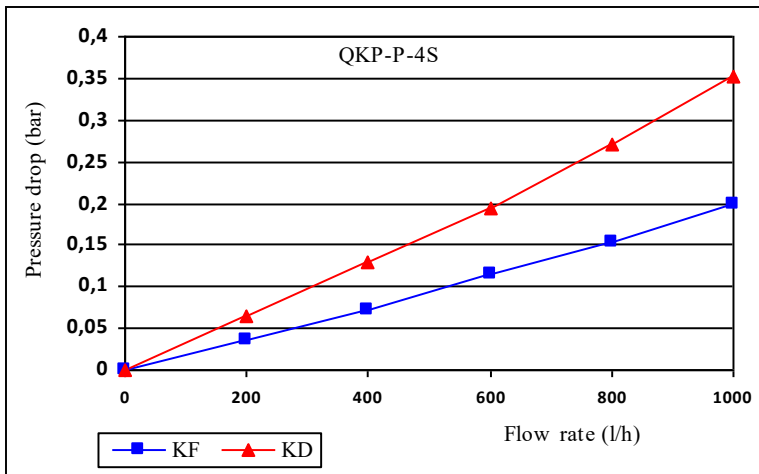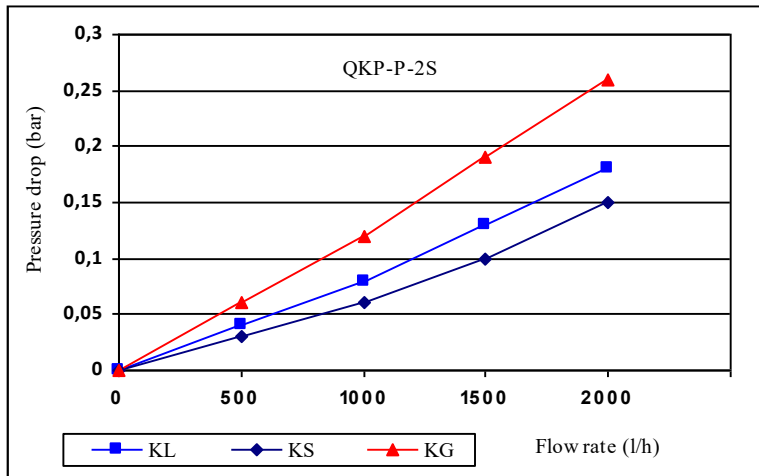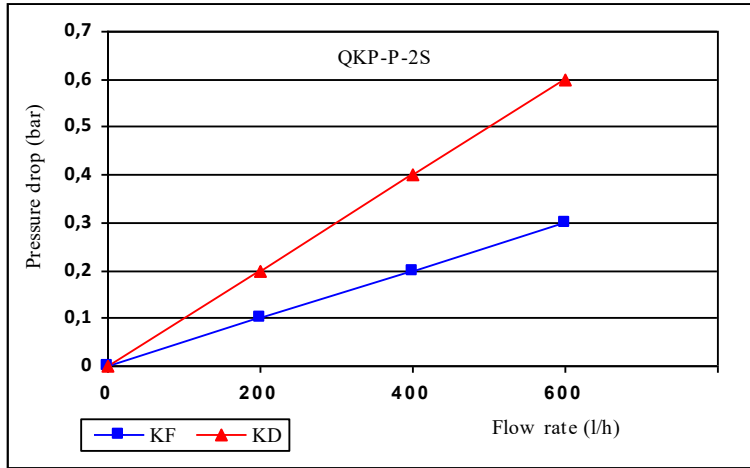

## QUALIKAP QKP-P ORDERING INFORMATION

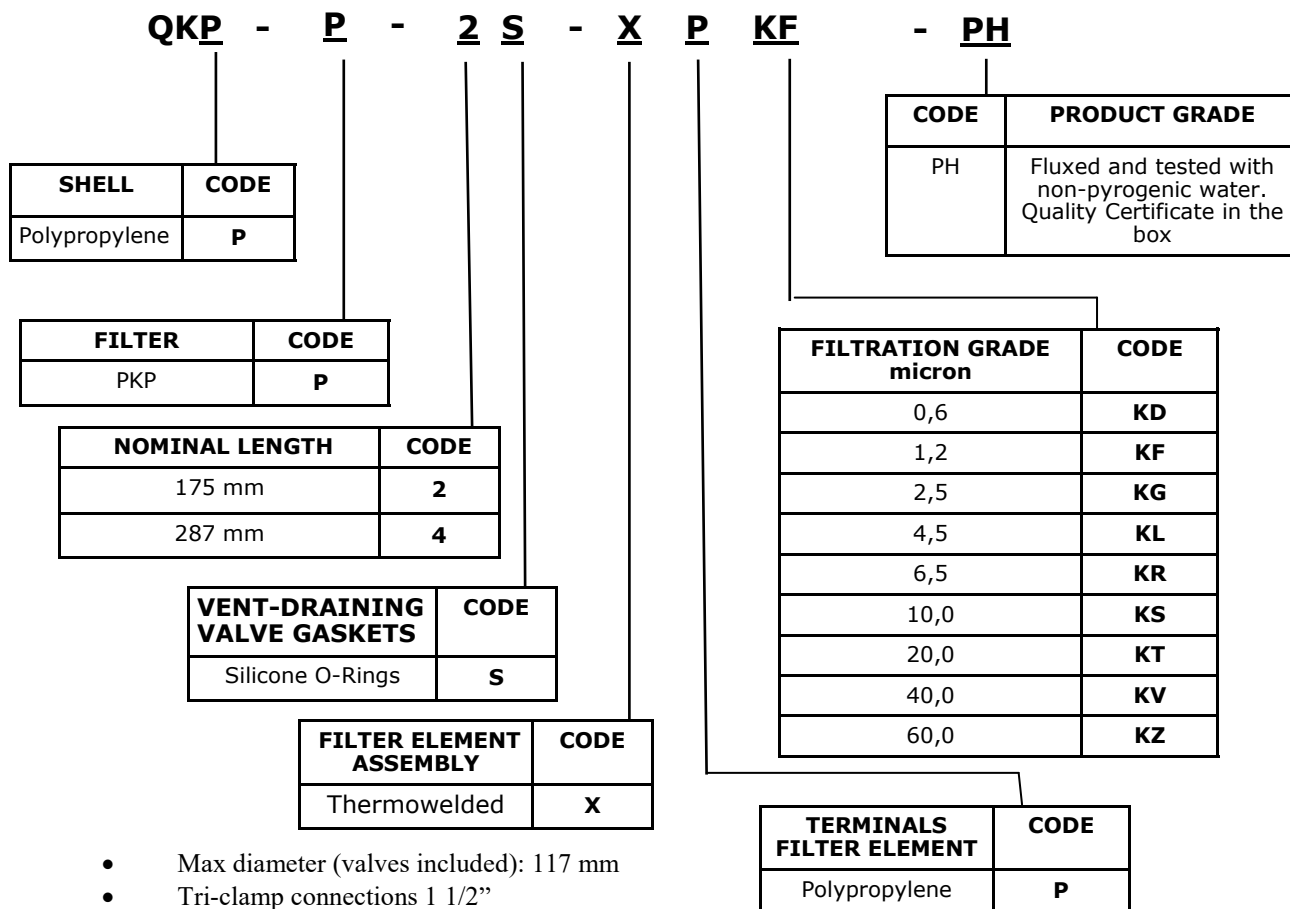

The data are informative and subject to change without notice. User is responsible for determining whether the product is fit for particular purpose and suitable for User's method of application.

### Bea Technologies Spa

Via Newton,4 -20016 Pero (MILANO) Italy  
 Tel.+(39) 02 339271 / Fax+(39) 02 3390713  
 mail:info@bea-italy.com  
 web:www.bea-italy.com

LY-QKPP-1019-EN-D6MG